



Nutri-Cath[®] with ENFit[™]

Nutri-Cath with ENFit is a family of locking enteral feeding-only catheters, extension sets, and oral dispensing syringes that:

- * Complies with ISO 80369-3
- * Utilizes the well-established GESCO[™] Nutri-Cath specialty enteral feeding catheters
- * Ensures the avoidance of a possible misconnection with devices administering medications
- * Ensures a secure connection that minimizes the risk of unwanted disconnections

30-Day Catheters

Nutri-Cath's gentle enteral-only catheter is a well-proven solution for either short term or intermittent feedings. Nutri-Cath resists irritation and prevents hyperactive gag reflex during nasogastric, orogastric or nasojejunal placement. Available in silicone and a thermosensitive hi-flow Tecoflex[®] polyurethane, the non-wettable materials will not perforate, clog or adhere to lining of the stomach.

Silicone Nutri-Cath with ENFit

Packaged 10/bx

	3.5 Fr	5.0 Fr	6.5 Fr	8.0 Fr
40 cm	4153517E	4155017E	4156517E	4158017E
58 cm	4153527E	4155027E	4156527E	4158027E
93 cm	4153537E	4155037E	4156537E	4158037E

Polyurethane Nutri-Cath with ENFit

Packaged 10/bx

	3.5 Fr	5.0 Fr	6.5 Fr	8.5Fr
40 cm	4153517PE	4155017PE	4156517PE	4158517PE
58 cm	4153527PE	4155027PE	4156527PE	4158527PE
93 cm	4153537PE	4155037PE	4156537PE	4158537PE

Extension Sets

Nutri-Cath with ENFit Continuous Feed Extension Sets

Packaged 50/bx

4150017E	60" Locking Extension Set, Female ENFit to Male ENFit
4150027E	60" Locking Extension Set w/ ENFit side port, Female ENFit to Male ENFit
4150067E	30" Locking Extension Set, Female ENFit to Male ENFit
415FIT	ENFit Transition Connector, Oral Female to Stepped Adaptor



Oral Dispensing Syringes



Monoject™ dispensing syringes:
For easy recognition by existing dispensing syringe pumps,
eliminating potential dosage errors

- * Purple plungers denote ENFit enteral feeding only
- * The secure locking ENFit connections are designed to reduce accidental disconnects and misconnects
- * Eliminates the need for a specialized enteral-only feeding pump

Sterile ENFit oral dispensing syringes, individually packaged

<u>Volume</u>	<u>Part #</u>	
1 mL:	41501E	60/bx
3 mL:	41503E	80/bx
6 mL:	41506E	80/bx
12 mL:	41512E	80/bx
35 mL:	41535E	40/bx
60 mL:	41560E	30/bx

Sterile ENFit oral dispensing syringe caps, 5 per package


	<u>Part #</u>	
CAP	415CPE	1000/bx

ISO 80369-3: 2016 specifies the dimensions and requirements for the design and functional performance of small-bore connectors intended to be used for connections on enteral medical devices and accessories. The ISO 80369-3 connector design in the reverse orientation is commonly known by its federally registered trademarked name ENFit.



UTAH MEDICAL PRODUCTS

FEMCARE

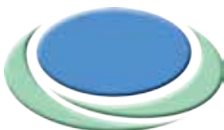

UNITED STATES

Utah Medical Products, Inc.
7043 South 300 West
Midvale, Utah 84047
Toll-free 800 533 4984
Tel (801) 566 1200
Fax (801) 566 2062
Email customer@utahmed.com
www.utahmed.com

AUSTRALIA

Femcare Australia
Unit 12, 5 Gladstone Road
Castle Hill NSW 2154
Tel +61 (02) 9045 4100
Fax +61 (02) 9894 1913
Email customer@femcare.com.au

EUROPE

Utah Medical Products Ltd.
Athlone Business & Technology Park
Garrycastle
Athlone, Co Westmeath N37 XK74
Ireland
Tel +353 (090) 647 3932
Fax +353 (090) 647 5608
Email customerservice_ireland@utahmed.com

Femcare Ltd.

32 Premier Way
Romsey, Hampshire S051 9DQ
United Kingdom
Tel +44 (01794) 525 100
Fax +44 (01794) 525 101
Email enquiries@femcare.co.uk